**bzmi**

SHM1 – toothbrush

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**Producent:** Shenzhen Minrek Electronic Ltd.


Adres: C602, 6th Floor, Building C, Hangcheng Road, Huafeng International Robot Industrial Zone, Xi'xiang, Bao'an, Shenzhen, Chiny 518126

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# USER MANUAL EN

Extended user guide available on [www.bemi-electronics.com](http://www.bemi-electronics.com)

## Product description

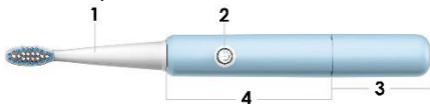1. Toothbrush head

2. Power button

3. Battery cover

4. Device main part

Ultrasonic Electric Toothbrush mode: SHM1, in manual called "device". It is mandatory to read this manual before start to using the device.

To **turn on the device**, press the function button. To turn off the device, press the power button for 2 seconds.

## Maintaining Equipment

Wash your brush head thoroughly before use. After each use, rinse the bristles and toothbrush head with clean water and air dry naturally. Clean the device regularly, keep it dry. Do not expose the device to cleaning agents, caustic agents or other chemicals. Do not expose the device to high (over + 35 ° C), low (below 0 ° C) temperatures or heat due to, for example, sunlight or other heat source. It may damage the device. The product is not suitable for children below 8 years and infants. Device can be used by children from 8 years old and above if they have been given supervision or instruction of concerning use of the device in a safe way, by an adult. **Note IPX7:** means that this device is resistant to dust, rain, moisture and a single immersion in standing water for a short time. Exposing the device to tougher conditions can damage it. When used for the first time, it may cause tooth discomfort, tingling and severe bleeding due to high-frequency vibration of the product. Please choose gentle mode, which will gradually disappear after repeated use. If there is


bleeding after two weeks of use, stop using this product and consult your dentist.

### **Battery change**

Note: Do not replace the battery when the device is wet. Do not allow moisture to get inside the device, for example when it is disassembled. The device is designed only to be powered by an AA battery, it is forbidden to use any other type of power source.

To replace the battery, follow the steps described in the photo at this point.

1. With one hand hold the main body of the device and with other hand turn the battery cover to the left, as in the photo below.
2. Slide off the battery cover
3. If some other battery is installed in the device, remove it by tilting it from the down part marked with the (+) sign. When the battery socket is empty, insert the new battery by first inserting the part marked with (-) sign. Make sure the battery fits.
4. Insert the battery cover into the device
5. Tighten the battery cover to the right as in the photo below.

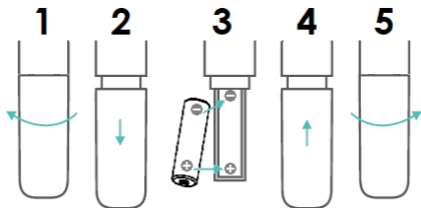
### **Toothbrush Head Replacement**


Its recommended to replace Toothbrush Head every 3 months or sooner if signs of wear appear.

When replacing hold the main part of the device with one hand and with the other hand rotate toothbrush head as on the picture bellow and pull in the opposite directions to the main

part of the device. Take out the toothbrush head, and then install a new toothbrush head. When installing new toothbrush head, please pay attention to the direction. The bristle direction of the toothbrush head should be in the same direction as the power button direction.

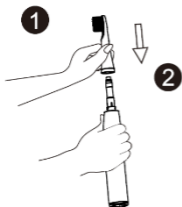ASSEMBLING

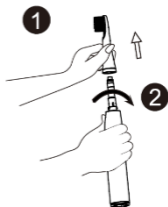DISASSEMBLING

**Basic Parameters:** Vibration frequency (per min.): 13 000, power source: replaceable AA battery; protection level: IPX7 waterproof.

**Proper disposal of this product.** The device cannot be disposed of together with household waste. To prevent potential threats to the environment and human health due to improper disposal of waste, it is mandatory to take used electronic equipment to a designated collection point, or to the place

where the product was purchased for the purpose of recycling electronic waste.

This device, Electric Toothbrush model SHM1, is in compliance with Directive 2014/30/EU of the European Parliament and of the Council. The complete declaration of conformity is available

online at the following link: [www.bemi-electronics.com/certificates](http://www.bemi-electronics.com/certificates)

This manual is not covered by the warranty, and the appearance and/or functionality of the software installed on the device may be different from the one described. The Guarantor reserves the right to modify the instruction manual at any time without prior notice.

**Manufacturer:** Shenzhen Minrek Electronic Ltd.

Address: C602, 6th Floor, Building C, Hangcheng Road, Huafeng International Robot Industrial Zone, Xi'xiang, Bao'an, Shenzhen, China 518126

Made in China