bzmi

SH1 – toothbrush

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Pobór mocy – 2W

Producent: Shenzhen Minrek Electronic Ltd.


Adres: C602, 6th Floor, Building C, Hangcheng Road, Huafeng International Robot Industrial Zone, Xi'xiang, Bao'an, Shenzhen, Chiny 518126

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USER MANUAL EN

Extended user guide available on www.bemi-electronics.com

Product description

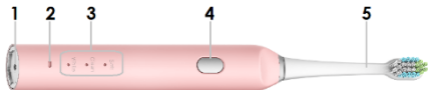1. Charging terminal 2. Power indicator 3. Cleaning modes
4. Function button 5. Toothbrush head

Ultrasonic Electric Toothbrush mode: SH1, in manual called "device". It is mandatory to read this manual before start to using the device.

To **turn on the device**, press the function button. To turn off the device, press the power button for 2 seconds.

Maintaining Equipment

Wash your brush head thoroughly before use. After each use, rinse the bristles and toothbrush head with clean water and air dry naturally. Clean the device regularly, keep it dry. Do not expose the device to cleaning agents, caustic agents or other chemicals. Do not expose the device to high (over + 35 ° C), low (below 0 ° C) temperatures or heat due to, for example, sunlight or other heat source. It may damage the device. The product is not suitable for children below 8 years and infants. Device can be used by children from 8 years old and above if they have been given supervision or instruction of concerning use of the device in a safe way, by an adult.

Note IPX7: means that this device is resistant to dust, rain, moisture and a single immersion in standing water for a short time. Exposing the device to tougher conditions can damage it.

Charging

Before start using the device, please make sure the battery is charged. For charging, please use the provided charging cable.

The recommended power adaptor current is 5V, 1A. It takes about 7 hours to fully charge device. Please follow the steps below to the charge.

1.Place the connector of the charging cable in the charging terminal of the device.

2.Connect the charging cable to the power source using the USB connector. The power indicator light (red light) start to flashing. Please charge the product for 4-7 hours or more in the following cases: Before the first use; Has not used this product for more than 6 months; In case of low power reminder.

When the power indicator is flashed red light, it indicates a low power. Insert a charging connection, with power, into device to charge when is turned on will automatically turn it off and will switch it in charging mode; During charging, the device cannot be turned on. When the device is fully charged, the power indicator will light green. After the device is used for a long time, the battery life will be shortened and the use time will be reduced.

Note: do not charge the device when it's wet, charging terminal, charging cable and power source must to be dry. Don't charge device close to watter.

Cleaning Modes

Programs can be changed by single press function button when device is turned on. Device has a memory function so when used, the product will remember the brushing mode of last shutdown.

1.Soft Mode: For sensitive teeth and gums.

2.Standard Mode: for routine cleaning and removing dental plaque.


3.Whitening Mode: Higher frequency, high speed tooth cleaning.

Smart timer (for every clining mode): The device has automatic cycles for easier control of cleaning time: pause every 30 seconds and after 2 minute turn-off.


When used for the first time, it may cause tooth discomfort, tingling and severe bleeding due to high-frequency vibration of

the product. Please choose gentle mode, which will gradually disappear after repeated use. If there is bleeding after two weeks of use, stop using this product and consult your dentist.

Toothbrush Head Replacement

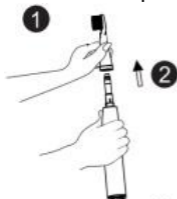Disassembly

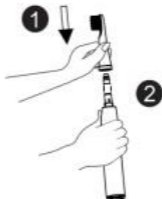installation

It is recommended to replace the Toothbrush Head every 3 months or sooner if signs of wear appear.

When replacing, hold the main part of the device tightly with one hand, hold the toothbrush head with the other hand, and pull it in opposite directions. Take out the toothbrush head, and then install a new toothbrush head. When replacing the toothbrush head, please pay attention to the direction. The bristle direction of the toothbrush head is in the same direction as the function button direction.

Basic Parameters: Vibration frequency (per min.): 36 000, battery capacity: 800 mAh; protection level: IPX7 waterproof.

Proper disposal of this product. The device cannot be disposed of together with household waste. To prevent potential threats to the environment and human health due to improper disposal of waste, it is mandatory to take used electronic equipment to a designated collection point, or to the place

where the product was purchased for the purpose of recycling electronic waste. **Warning:** Lithium battery inside.

This device, Electric Toothbrush model SH1, is in compliance with Directive 2014/30/EU of the European Parliament and of the Council. The complete declaration of conformity is available

online at the following link: www.bemi-electronics.com/certificates

This manual is not covered by the warranty, and the appearance and/or functionality of the software installed on the device may be different from the one described. The Guarantor reserves the right to modify the instruction manual at any time without prior notice.

Power consumption: 2W

Manufacturer: Shenzhen Minrek Electronic Ltd.

Address: C602, 6th Floor, Building C, Hangcheng Road, Huafeng International Robot Industrial Zone, Xi'xiang, Bao'an, Shenzhen, China 518126

Made in China